

Midi Plus™

Эргономичное дозирование больших объемов



Легкое и эргономичное дозирование больших объемов из бутылей или пробирок


Эргономичный дизайн позволит избежать RSI

Midi Plus™ - это электронное устройство управления дозированием (насос), небольшого веса, позволяющее набирать жидкости из бутылей и пробирок. Эргономичная конструкция устройства помогает свести к минимуму нагрузку на руку оператора, которая неизбежна в случае использования серологических или градуированных пипеток. Благодаря электронному управлению работой Midi Plus™, усилие пальца, необходимое для наполнения дозатора, гораздо меньше, чем у устройств, использующих сжимающуюся грушу. Эти особенности Midi Plus™ позволяют предотвратить утомление или травмирование руки оператора, а также снизить риск развития синдрома лучезапястного канала (RSI - Repetitive Strain Injury).



Плавная регулировка скорости обеспечивает точность дозирования

- Скорость набора и дозирования жидкости легко регулировать с помощью регулировочного колеса или нажатия на операционные кнопки.
- Выбор скорости работы насоса может осуществляться в зависимости от особенностей применения и объема пипетки.
- Возможность гравитационного (капельного) дозирования. Просто благодаря меньшему усилию при нажатии кнопки дозирования.

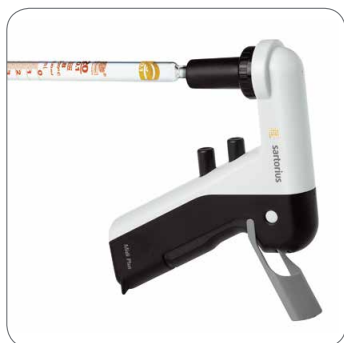

Экономия места за счет откидной подставки

- Встроенная откидная подставка позволяет сэкономить место на рабочем столе.
- Когда Midi Plus™ стоит на подставке, стеклянная или пластиковая пипетка не соприкасается с поверхностью стола, что предотвращает возможную контаминацию.





Совместимость с большинством пипеток и наконечников

- Насос Midi Plus™ совместим с большинством стеклянных и пластиковых пипеток объемом 1 – 100 мл.
- Насос Midi Plus™ также совместим со сменными наконечниками Sartorius Optifit объемом 5 мл и 10 мл, которые можно использовать, когда точность дозирования не существенна.



Гидрофобные фильтры предотвращают набор избыточного объема жидкости

- Гидрофобные фильтры Sartorius предотвращают набор избыточного объема жидкости, благодаря чему исключается риск контаминации насоса.
- Фильтры поставляются как в стерильном, так и в нестерильном исполнении, в зависимости от требования режима чистоты.

Информация для оформления заказа. Аксессуары.

Номер по каталогу(код)	Наименование	Кол-во/ Ед.Иtem
710931	Midi Plus™	1
LN-7129120	Сменный фильтр 0,45 мкм, нестерильный, автоклавируемый	5
LN-7129130	Сменный фильтр, 0,2 мкм, предварительно стерилизованный	5
LN-711017	Адаптер для наконечников 5 мл, автоклавируемый, включая конусный наконечник и силиконовый адаптер	1
LN-711018	Адаптер для наконечников 10 мл, автоклавируемый, включая конусный наконечник и силиконовый адаптер	1
780300	Наконечник Optifit 5 мл	100
780308	Наконечник Optifit 5 мл	1000
780310	Наконечник Optifit 10 мл (удлинённый для Midi Plus™)	100



Sales and Service Contacts

For further contacts, visit www.sartorius.com

Europe

Germany

Sartorius Lab Instruments
GmbH & Co. KG
Weender Landstrasse 94-108
37075 Goettingen
Phone +49.551.308.0
Fax +49.551.308.3289

France & Suisse Romande

Sartorius France S.A.S.
ZI Les Paluds
Avenue de Jouques - CS 31090
13781 Aubagne Cedex
Phone +33.442.845600
Fax +33.442.846545

Austria

Sartorius Austria GmbH
Franzosengraben 12
1030 Vienna
Phone +43.1.7965760.0
Fax +43.1.7965760.24

Belgium

Sartorius Belgium N.V.
Leuvensesteenweg, 248/B
1800 Vilvoorde
Phone +32.2.756.06.71
Fax +32.2.253.45.95

Finland & Baltics

Sartorius Biohit Liquid Handling Oy
Laippatie 1
00880 Helsinki
Phone +358.9.755.951
Fax +358.9.755.95.292

Hungary

Sartorius Hungária Kft.
Kagyló u. 5.
2092 Budakeszi
Phone +3623.457.227
Fax +3623.457.147

Ireland

Sartorius Ireland Ltd.
Unit 41, The Business Centre
Stadium Business Park
Ballycoolin Road
Dublin 11
Phone +353.1.8089050
Fax +353.1.8089388

Italy

Sartorius Italy S.r.l.
Via dell' Antella, 76/A
50012 Antella-Bagno a Ripoli (FI)
Phone +39.055.63.40.41
Fax +39.055.63.40.526

Poland

Sartorius Poland sp.z o.o.
ul. Wrzesinska 70
62-025 Kostrzyn
Phone +48.61.6473830
Fax +48.61.6473839

Russian Federation

LLC "Sartorius ICR"
and LLC "Biohit"
Uralskaya str. 4, Lit. B
199155, Saint-Petersburg
Phone +7.812.327.5.327
Fax +7.812.327.5.323

Spain & Portugal

Sartorius Spain S.A.
Offices in Madrid:
C/ Isabel Colbrand, 10-12, of. 70
28050 Madrid
Phone Spain +34.902.123.367
Phone Portugal +351.800.855.800
Fax Spain +34.91.358.84.85
Fax Portugal +351.800.855.799

Switzerland

Sartorius Mechatronics Switzerland AG
Ringstrasse 24a
8317 Tagelswangen (ZH)
Phone +41.44.746.50.00
Fax +41.44.746.50.50

U.K.

Sartorius UK Ltd.
Longmead Business Centre
Blenheim Road, Epsom
Surrey KT19 9QQ
Phone +44.1372.737159
Fax +44.1372.729927

America

USA

Sartorius Corporation
5 Orville Drive, Suite 200
Bohemia, NY 11716
Phone +1.631.254.4249
Toll-free +1.800.635.2906
Fax +1.631.254.4253

Argentina

Sartorius Argentina S.A.
Int. A. Avalos 4251
B1605ECS Munro
Buenos Aires
Phone +54.11.4721.0505
Fax +54.11.4762.2333

Brazil

Sartorius do Brasil Ltda
Avenida Senador Vergueiro 2962
São Bernardo do Campo
CEP 09600-000 - SP- Brasil
Phone +55.11.4362.8900
Fax +55.11.4362.8901

Canada

Sartorius Canada Inc.
2179 Dunwin Drive #4
Mississauga, ON L5L 1X2
Phone +1.905.569.7977
Toll-Free +1.800.668.4234
Fax +1.905.569.7021

Mexico

Sartorius de México S.A. de C.V.
Circuito Circunvalación Poniente
No. 149
Ciudad Satélite
53100, Estado de México
México
Phone +52.5555.62.1102
Fax +52.5555.62.2942

Asia|Pacific

Australia

Sartorius Australia Pty. Ltd.
Unit 5, 7-11 Rodeo Drive
Dandenong South Vic 3175
Phone +61.3.8762.1800
Fax +61.3.8762.1828

China

Sartorius Scientific
Instruments (Beijing) Co., Ltd.
33 Yu An Road, Airport Industrial Park
Zone B,
Shunyi District, Beijing 101300,
P.R.China
Phone +86.10.8042.6300
Fax +86.10.8042.6486

Hong Kong

Sartorius Hong Kong Ltd.
Unit 1012, Lu Plaza
2 Wing Yip Street
Kwong Tong
Kowloon, Hong Kong
Phone +852.2774.2678
Fax +852.2766.3526

India

Sartorius Weighing India Pvt. Ltd.
#69/2-69/3, NH 48, Jakkasandra,
Nelamangala Tq
562 123 Bangalore, India
Phone +91.80.4350.5250
Fax +91.80.4350.5253

Japan

Sartorius Stedim Japan K.K.
4th Fl., Daiwa Shinagawa North Bldg.
8-11, Kita-Shinagawa 1-chome
Shinagawa-ku, Tokyo, 140-0001 Japan
Phone +81.3.4331.4300
Fax +81.3.4331.4301

Malaysia

Sartorius Malaysia Sdn. Bhd
Lot L3-E-3B, Enterprise 4
Technology Park Malaysia
Bukit Jalil
57000 Kuala Lumpur, Malaysia
Phone +60.3.8996.0622
Fax +60.3.8996.0755

Singapore

Sartorius Singapore Pte. Ltd
1 Science Park Road,
The Capricorn, #05-08A,
Singapore Science Park II
Singapore 117528
Phone +65.6872.3966
Fax +65.6778.2494

South Korea

Sartorius Korea Ltd.
8th Floor, Solid Space B/D,
PanGyoYeok-Ro 220, BunDang-Gu
SeongNam-Si, GyeongGi-Do
463-400
Phone +82.31.622.5700
Fax +82.31.622.5799

Thailand

Sartorius (Thailand) Co. Ltd.
129 Rama 9 Road,
Huaykwang
Bangkok 10310
Phone +66.2643.8361-6
Fax +66.2643.8367



◀ www.sartorius.com